

Comfort-Sonic®

LED Piezo Ultrasonic Scaler

Simply Practical

Boost your productivity and enhance patient experience with Comfort-Sonic®, the LED Piezo Scaler that is as versatile as your operational needs. It is compact and easy to use for a variety of daily applications including: scaling, endo and perio. Effortlessly eliminate calculus and other stubborn buildup with its ergonomic lightweight handpiece that features six LED lights encircling the tip for optimal illumination. The precise power delivered through steady linear ultrasonic vibrations ensures not only low operating noise, but also that movement is concentrated only at the tip, resulting in handpiece stability and maximized user comfort.



Powerful, with easy adjustable design



Built-in water regulation



3 Modes: scaling, perio, endo



Detachable handpiece



Autoclavable up to 135°C




Daylight LED: maximize visibility & minimize eye fatigue

Key Specifications

Main unit input	24V~ 50Hz/60Hz 1.3A
Output power	3W ~ 20W
Frequency	28kHz ± 3kHz
Power input	110V~ 50Hz/60Hz, 220 - 240V~ 50Hz/60Hz
Water pressure	0.1bar ~ 5bar(0.01MPa ~ 0.5MPa)
Main unit weight	P6 Plus: 0.64kg, P9:1.35kg
Dimensions	P6 Plus: 196mm×134mm×80mm, P9: 234mm×193mm×109mm

Order Information

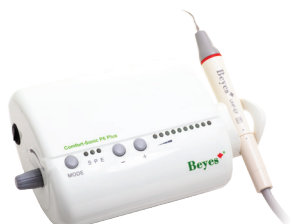

Connect directly to your waterline

Model: **P6 Plus**

Code: **UL2104**

Intro kit includes:

- 6 Tips (2 EG1, EG2, EG4, EP1, EE1)
- 4 NITI U-Files
- 1 Handpiece
- 1 WST Wrench for scaler tips
- 1 WEF Wrench for endo files
- 1 Power supply
- 1 Foot pedal



Optimized portability

Model: **P9**

Code: **UL2103**

Intro kit includes:

- 8 Tips (EG1, EG4, EG5, EG6, EP1, EP3, EP4D, EE1)
- 2 Solution containers
- 1 Tip holder
- 4 NITI U-Files
- 1 Handpiece
- 1 WST Wrench for scaler tips
- 1 WEF Wrench for endo files
- 1 Power supply
- 1 Foot pedal



Accessories and Tips

Option 1: Scan the bar code above for the Tips & Accessories brochure.

Option 2: Go to the website link:
<http://beyes.ca/en/categories/19/products/43>



Beyes Dental Canada Inc.
23-595 Middlefield Road
Toronto, Ontario, M1B 3M5
Canada
Tel: 1-855-603-1888
Email: info@beyes.ca
Web: www.beyes.ca



Lotus NL B.V.
Koningin Julianaplein 10, le Verd,
2595AA, The Hague,
Netherlands
Tel: +31645171879 (English)
+31626669008 (Dutch)

Your Dealer for Beyes Products

Scan the QR code for more information on Comfort-Sonic



Printed in Canada
File #.ENB015R2

© Beyes Canada Inc 2019. All Rights Reserved.
Specifications are subject to change without notice.
All brands are holders of their respective trademarks.
All Beyes products are manufactured in strict compliance with Health Canada regulations ensuring quality.
Manufacturing facility is ISO 13485 certified under Canadian Medical Devices Conformity Assessment System (CMDCAS).