

# AirLight<sup>®</sup> M800 Smart

*No Light Bulb*

*No Fiber Optic Tubing*

*No Electrical Source*

*But you will get **LED Illumination***

The newest AirLight will exceed any previous market expectations on LED illumination. The truly innovative design of this handpiece incorporates a highly efficient micro generator to power its LED. It will eliminate the need of fiber optic system from your dental delivery unit. This cutting edge technology allows the operators to have a clearer illuminated operative field to reduce fatigue. The Direct-LED is brighter and has larger light pattern than the traditional fiber optic system. The well designed and precisely manufactured AirLight provides extraordinary torque to ensure the powerful cutting that you need. It will fit directly into your existing 4, 5, 6 holes handpiece tubings.



**Award  
Winning:**



## Key Features



### Direct-LED

LED is optimally positioned on the face of the handpiece to provide a larger light pattern than traditional fiber optic systems.



### Generator

An innovative micro generator that produces electricity for the LED is built into the handpiece. As a result, it eliminates the need for a fiber optic system.



### Precision Ball Bearings

It is specially designed for smooth and quiet operation. It also provides superior torque and a long lifespan.




### Micro-Tex Surface

The latest surface enhancement technology creates a microtextured and scratch proof surface that lowers the possibility of slipping during operation.



### Triple Water Spray

Effectively cools the entire surface of the bur that avoids overheating and provides better patient comfort.

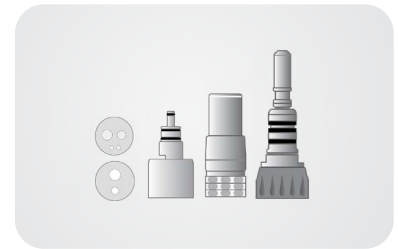

### Compatibility

Compatible with a great variety of coupling systems and adapts well to every existing treatment unit.

## Key Specification

Colour temperature (K)	5,500
Brightness (LUX)	25,500
Head Size (Width x Height)	10.8 x 13.0mm
Drive Air Pressure(PSI)	40-60
Power (WATT)	28
Speed (RPM)	410,000
Rotary instruments	FG burs according to EN ISO 1797-1
Ø / max. length (MM)	1.6 / 25
Water Spray	Triple Triple



## Advantages

AirLight offers many advantages compared to the traditional fiber optic whether it has halogen or LED bulb. The light bulb of traditional fiber optic located on the tubing or the coupler. The light is then conducted by a fiber optic rod to the head of the handpiece. As a result, the intensity is reduced during the transition and the coverage is limited. Additionally, the fiber optic rod gets degraded over time with autoclaving, and your light becomes getting dimmer and yellowish. Eventually, you will have to replace the fiber optic rod.



Fiber Optic  
Halogen Bulb

Fiber Optic  
LED Bulb

**AirLight®**  
Direct-LED

**AirLight®**  
**SAVES YOUR MONEY**  
on:

- Fiber optic systems
- Replacement tubings
- Replacement bulbs
- Maintenance costs

*AirLight is completely different and has none of the above problems because the LED is located on the face of the handpiece to produce direct illumination. The electricity is supplied by two medical graded cables from the generator. The LED is engineered in Germany with macro angle design and high colour rendering index. As a result, you will get higher intensity and a larger coverage.*

*"...It does not just illuminate a small spot, it is the entire oral cavity..." -Dr. Gaw*

## Order Information

### 4 Hole Backend



**Code:** HP3008  
**Model:** M800-M/M4

### 2 Hole Backend



**Code:** HP3051  
**Model:** M800-M/B2

### Beyes® PD Backend



**Code:** HP6018  
**Model:** M800-M/PD

### NSK® Backend



**Code:** HP3043  
**Model:** M800-M/N

### Kavo® Backend



**Code:** HP3036  
**Model:** M800-M/K

### EZ-Swivel® Coupler



**Code**    **Model**  
HP6008    PD M4  
HP6002    PD B2



**Code**    **Model**  
HP2063    K M4  
HP2064    K E6

According to thousands of users,  
AirLight is proven to be the future of Dental handpiece!



*"I feel lucky to have found and be able to work with the new AirLight V3 from Beyes. This feels great in my hand: the V3 has a new smart feature a bright self-powered LED light. I found it brighter than my other actual LED fiber optic handpieces. I am happy to contribute to the success of a Canadian company. For all these reasons I recommend the AirLight V3 without any hesitations."*

**Dr. Matthieu Menard**  
QC, Canada



*"It's unique and very bright feature of led light without the need for fiber optic tubing seemed to me attractive. The powerful water spray and bright illumination coupled with high torque value of the hand pieces and the very convenient swivel coupler all help out in our daily routine of practicing quality dentistry in a comfortable environment. I found this to be very user friendly and forthcoming."*

**Dr. Joshua Haimovici**  
QC, Canada



*"We have ordered many handpieces from Beyes Canada Inc. The handpieces are light, smooth, and reliable. The LED handpiece increases visibility, and decreases strain on the eyes which causes fatigue. The handpieces are definitely a good value, and we highly recommend that you integrate them into your office as well."*

**Dr. Ann Sahota**  
BC, Canada



*"The incorporation of the Beyes AirLight high speed handpiece has proven a valuable tool for our clinical armamentarium, with great advantages of the light weight, high torque and brightness of the fiber-optics compared to other handpieces. The hand pieces work great, have good power and water spray. Their products are great and their prices are way better than other manufacturers"*

**Dr. Jaspal Singh Sarao**  
B.C, Canada



*"The Air Light M800-S, despite the amount of usage and frequency of autoclaving, never seems to have a compromise in the department of torque, durability, and dependability. Each and every cut is just as precise as the one before."*

**Dr. Frank M. Glusheski**  
PA, USA



*"I was impressed by the new technology of your hand-piece, a fiber-optic handpiece without an external light source and without the cumbersome fiber-optic hose. I can say that this is one of the best dental purchases that I have ever made. This handpiece is competitive, if not better than any of the best branded handpieces. For the money, it is the best value out in the market."*

**Dr. Allen Gaw**  
TX, USA



*"The Airlight's light, fits well in hand and easy to operate. The water spray is awesome, cooling down effectively the bur and at the same time does not fill the patient mouth with water. The fibreoptic system is even better, light is bright and white giving a wider vision and not focused on a single spot. Having all these cool features in a single handpiece is really great."*

**Dr. Pierre Hayek**  
QC, Canada



*"I love our AirLight handpieces. We've had no problems since we bought them and they have great torque and water flow. It's swivelling connector makes it comfortable to hold. The light is especially beneficial in those hard to see areas, even when I use my headlamp. The AirLight handpiece offers tremendous value along with great performance."*

**Dr. Neel Joshi**  
ON, Canada



*"It was very light and easy to maneuver."*

**Dr. Frank Kuzmin**  
CT, USA



*"Thanks for introducing me to these amazing handpieces! They are now in all of our offices!"*

**Dr. Charles Botbol**  
ON, Canada

**Beyes**  
The Next Standard

Beyes Dental Inc.  
55 Laurier Drive  
Morrisburg, Ontario  
K0C 1X0, Canada  
Tel: 1-855-603-1888  
Email: info@beyes.ca  
Website: www.beyes.ca



**Lotus NL B.V.**  
Koningin Julianaplein 10, le Verd,  
2595AA, The Hague,  
Netherlands  
Email: info@lotusnl.com

Your Dealer for Beyes Products



Scan the QR code  
for more information

Printed in Canada  
File #: ENB012R2  
Rev.3/11.12.25

© Beyes Canada Inc 2025. All Rights Reserved.  
Specifications are subject to change without notice.  
All brands are holders of their respective trademarks.  
All Beyes products are manufactured in strict compliance  
with Health Canada regulations ensuring quality.  
Manufacturing facility is ISO 13485 certified under the  
Medical Device Single Audit Program (MDSAP).