


AirLight®
M800

X

traordinary

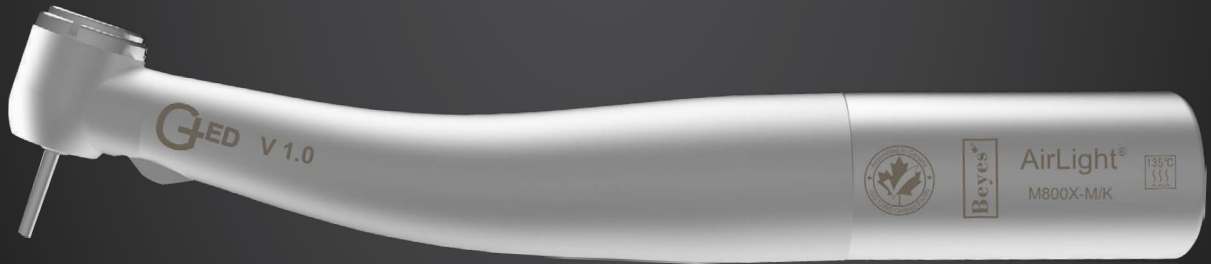

POWERFUL **Q**UIET **S**AFE **C**OMFY

Key Features

X Ball Bearings

Newly engineered ball bearings with special caging material that lasts four times longer than other ball bearings.

Limitless PSI Technology

Experience the freedom from the need for operating pressure control

Quattro Spray

Four specially-positioned water ports guaranteeing optimal cooling.

35W Torque

Increases efficiency and minimizes stall for the best precision cutting.

Direct-LED+

Provides brighter and larger light pattern enabling clinicians to have a clearer operative field.

IS Technology

Provides instant stop capability which eliminates backflow and noise when ceasing operation.

Micro-Tex

The slip-proof coating that improves hand grip, adding safety and comfort.

Built-in Generator

Eliminates the need for being tethered to a delivery unit's fiber optic system.

The well-designed and precisely-manufactured AirLight M800X delivers an extraordinary torque of (35W) ensuring accurate cutting. It is equipped with a specially designed (X-Ball Bearings) for quieter sound and longer life. The handpiece also possesses instant-stop capability (IS Tech) that lessens noise and eliminates backflow. In addition to four specially-positioned water ports guaranteeing optimal cooling, it has an all-new (Direct-LED+) that provide brighter and larger light pattern enabling clinicians to have a clearer operative field. Finally, the new handpiece comes with an extended (25 Months) Warranty.



M800X

Beyes® Handpieces

Model	AirLight® M800X	AirLight® M800S	AirLight® M800P	Maxso® X200S	Maxso® M200S	Maxso® M200E
RPM	430,000	400,000	400,000	400,000	400,000	400,000
Power	35	28	28	28	23	20
Spray Ports	4	3	3	3	3	1
Illumination	LED+	LED	LED	Fiber Optic	—	—
Instant Stop	✓	—	—	—	—	—
Built-in Generator	✓	✓	✓	—	—	—
Material	Stainless Steel	Special Alloy	Titanium	Special Alloy	Special Alloy	Special Alloy
Warrantee	25	12	12	12	12	6

Direct-LED



VS

Fiber Optic



Award Winning AirLight® Family



M800X

Key Specification

Power (WATT)	35
Drive Air Pressure (PSI)	40-60
Speed (RPM)	430,000 (10%±)
Water Spray	Quattro
Colour Temperature (K)	5,500
Brightness (LUX)	25,500
Ø / max. length (MM)	1.6 / 25
Chuck Type	Push Button
Head Size (W x H) mm	10.8 x 13.0
Rotary Instruments	FG burs according to EN ISO 1797-1

Order Information

Beyes® Backend



Code: **HP6018X**
Model: **M800X-M/PD**

NSK® Backend



Code: **HP3043X**
Model: **M800X-M/N**

Kavo® Backend



Code: **HP3036X**
Model: **M800X-M/K**

Star® Backend



Code: **HP3054X**
Model: **M800X-M/ST**



Beyes Dental Inc.
55 Laurier Drive
Morrisburg, Ontario
K0C 1X0, Canada
Tel: 1-855-603-1888
Email: info@beyes.ca
Website: www.beyes.ca



Lotus NL B.V.
Koningin Julianaplein 10, le Verd,
2595AA, The Hague,
Netherlands
Email: info@lotusnl.com

Your Dealer for Beyes Products



Scan the QR code
for more information

Printed in Canada
File #: ENB004
Rev.7/11.12.25

© Beyes Canada Inc 2025. All Rights Reserved.
Specifications are subject to change without notice.
All brands are holders of their respective trademarks.
All Beyes products are manufactured in strict compliance
with Health Canada regulations ensuring quality.
Manufacturing facility is ISO 13485 certified under the
Medical Device Single Audit Program (MDSAP).