


AirLight®
M800

X

traordinary

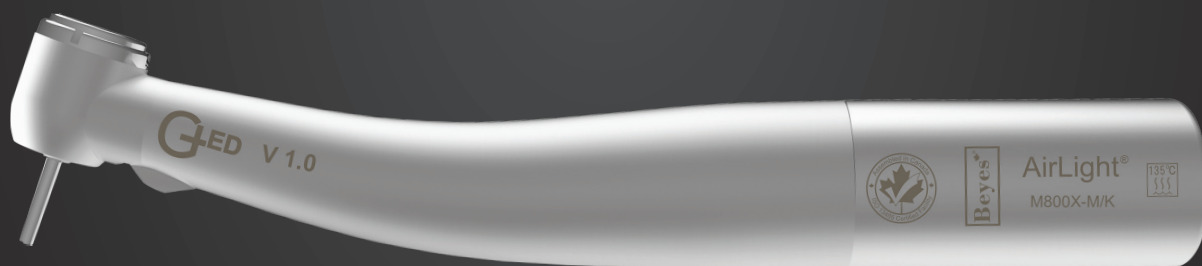

POWERFUL QUIET SAFE COMFY

Key Features

27W Torque

Increases efficiency and minimizes stall for the best precision cutting.

X Ball Bearings

Newly engineered ball bearings with special caging material that lasts four times longer than other ball bearings.

Quattro Spray

Four specially-positioned water ports guaranteeing optimal cooling.

IS Technology

Provides instant stop capability which eliminates backflow and noise when ceasing operation.

Direct-LED+

Provides brighter and larger light pattern enabling clinicians to have a clearer operative field.

Built-in Generator

Eliminates the need for being tethered to a delivery unit's fiber optic system.

Micro-Tex

The slip-proof coating that improves hand grip, adding safety and comfort.

The well-designed and precisely-manufactured AirLight M800X delivers an extraordinary torque of (27W) ensuring accurate cutting. It is equipped with a specially designed (X-Ball Bearings) for quieter sound and longer life. The handpiece also possesses instant-stop capability (IS Tech) that lessens noise and eliminates backflow. In addition to four specially-positioned water ports guaranteeing optimal cooling, it has an all-new (Direct-LED+) that provide brighter and larger light pattern enabling clinicians to have a clearer operative field. Finally, the new handpiece comes with an extended (25 Months) Warranty.



M800X

Beyes® Handpieces

Model	AirLight® M800X	AirLight® M800S	AirLight® M800P	Maxso® X200S	Maxso® M200S	Maxso® M200E
RPM	430,000	400,000	400,000	400,000	400,000	400,000
Power	27	23	23	23	21	21
Spray Ports	4	3	3	3	3	1
Illumination	LED+	LED	LED	Fiber Optic	—	—
Instant Stop	✓	—	—	—	—	—
Built-in Generator	✓	✓	✓	—	—	—
Material	Stainless Steel	Special Alloy	Titanium	Special Alloy	Special Alloy	Special Alloy
Warranty	25	12	12	12	12	6

Direct-LED



VS

Fiber Optic



Award Winning AirLight® Family



M800X

Key Specification

Power (WATT)	27
Drive Air Pressure (PSI)	40-60
Speed (RPM)	430,000 (10%±)
Water Spray	Quattro
Colour Temperature (K)	5,500
Brightness (LUX)	25,500
Ø / max. length (MM)	1.6 / 25
Chuck Type	Push Button
Head Size (W x H) mm	10.8 x 13.0
Rotary Instruments	FG burs according to EN ISO 1797-1

Order Information

Beyes® Backend



Code: **HP6018X**
Model: **M800X-M/PD**



NSK® Backend



Code: **HP3043X**
Model: **M800X-M/N**

Kavo® Backend



Code: **HP3036X**
Model: **M800X-M/K**

Star® Backend



Code: **HP3054X**
Model: **M800X-M/ST**

Beyes®
The Next Standard

Beyes Dental Canada Inc.
23-595 Middlefield Road
Toronto, Ontario, M1V 3S2
Canada
Tel: 1-416-281-1888
Email: info@beyes.ca
Web: www.beyes.ca

Your Dealer for Beyes Products

Scan the QR code for
more information on M800X



Printed in Canada
Document #: ENB004
Rev.3 / 06.08.19

© Beyes Canada Inc 2019. All Rights Reserved.
Specifications are subject to change without notice.
All brands are holders of their respective trademarks.
All Beyes products are manufactured in strict compliance
with Health Canada regulations ensuring quality.
Manufacturing facility is ISO 13485 certified under Canadian
Medical Devices Conformity Assessment System (CMDACS).