


Fiber-optic

Maxso[®] X200 Smart

The Maxso series high speed handpiece is designed to meet expectations in modern dental procedures. High quality components and comfortable grip with Micro-tex finish provide superior handling that you desire. Furthermore, the well aimed triple water spray will provide excellent coolant effects.

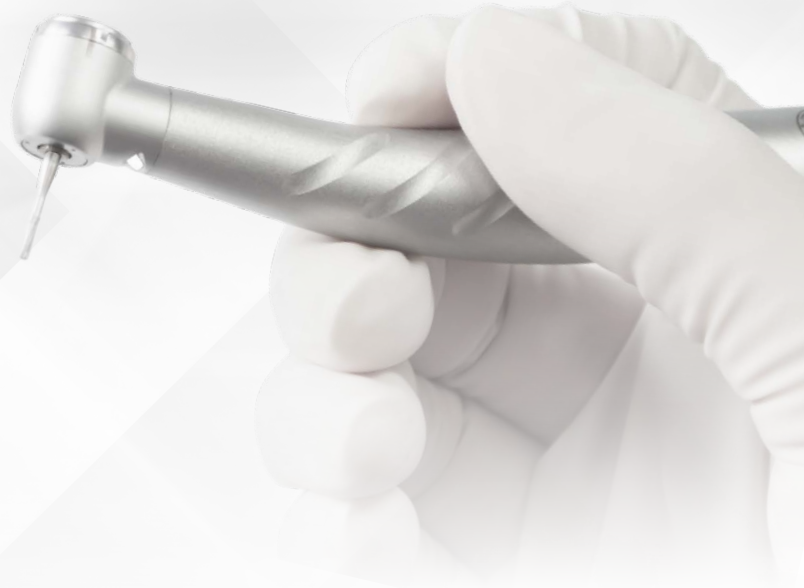


Triple Spray

Effectively cools the entire surface of the bur that avoids overheating and provides better patient comfort.



Ceramic

It is specially designed for smooth and quiet operation. It also provides superior torque and a long lifespan.

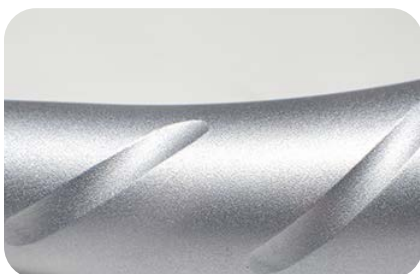

Micro-Tex

The latest surface enhancement technology creates a micro-textured and scratch proof surface that lowers the possibility of slipping during operation.



Power

Strong chuck holding force to prevent slippage and improve power with 23W torque.

Key Specification

	Mini Head	Standard Head
Head Size (Width x Height)	10.8mm x 13.0mm	12.5mm x 13.1mm
Drive Air Pressure (PSI)	40-60 PSI	40-60 PSI
Power (WATT)	21	23
Speed (RPM)	400,000	350,000
Rotary Instruments	FG burs according to EN ISO 1797-1	
Ø / max. length (MM)	1.6 / 25	1.6 / 25
Chuck Type	Push Button	Push Button
Water Spray	Triple	Triple

Order Information

Kavo® Backend



Mini Head

Code: **HP2060**
Model: **X200 M/K**

Standard Head

Code: **HP2061**
Model: **X200 S/K**

NSK® Backend



Mini Head

Code: **HP2066**
Model: **X200 M/N**

Standard Head

Code: **HP2067**
Model: **X200 S/N**



Beyes
CANADA

Beyes Dental Canada Inc.
23-595 Middlefield Rd
Scarborough, ON M1V 3S2
Canada
Tel: 1-855-603-1888
Email: info@beyes.ca
Website: www.beyes.ca



Lotus NL B.V.
Koningin Julianaplein 10, le Verd,
2595AA, The Hague,
Netherlands

Your Dealer for Beyes Products



Scan the QR code
for more information

Printed in Canada
File #: ENB011R2
Rev.2/09.16.24

© Beyes Canada Inc 2024. All Rights Reserved.
Specifications are subject to change without notice.
All brands are holders of their respective trademarks.
All Beyes products are manufactured in strict compliance
with Health Canada regulations ensuring quality.
Manufacturing facility is ISO 13485 certified under the
Medical Device Single Audit Program (MDSAP).