


Fiber-optic

Maxso®


X200 Smart

The Maxso series high speed handpiece is designed to meet expectations in modern dental procedures. High quality components and comfortable grip with Micro-tex finish provide superior handling that you desire. Furthermore, the well aimed triple water spray will provide excellent coolant effects.

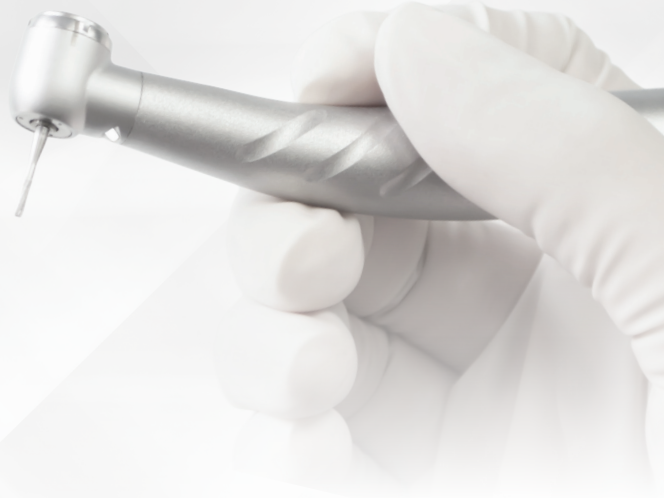


Triple Spray

Effectively cools the entire surface of the bur that avoids overheating and provides better patient comfort.

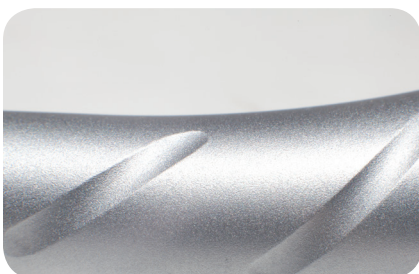

Micro-Tex

The latest surface enhancement technology creates a micro-textured and scratch proof surface that lowers the possibility of slipping during operation.

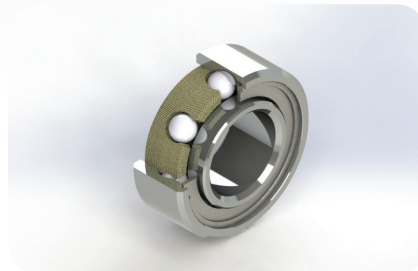

Ceramic

It is specially designed for smooth and quiet operation. It also provides superior torque and a long lifespan.



Power

Strong chuck holding force to prevent slippage and improve power with 23W torque.

Technical Data

	Mini Head	Standard Head
Head Size (Width x Height)	10.8mm x 13.0mm	12.5mm x 13.1mm
Drive Air Pressure (PSI)	40-60 PSI	40-60 PSI
Power (WATT)	21	23
Speed (RPM)	400,000	350,000
Rotary Instruments	FG burs according to EN ISO 1797-1:1995	
Ø / max. length (MM)	1.6 / 25	1.6 / 25
Chuck Type	Push Button	Push Button
Water Spray	Triple	Triple

Order Information

Kavo® Connection

**Mini Head**

Code **HP2060**
Model **X200 M/K**

Standard Head

Code **HP2061**
Model **X200 S/K**

NSK® Connection

**Mini Head**

Code **HP2066**
Model **X200 M/N**

Standard Head

Code **HP2067**
Model **X200 S/N**



Beyes®

The Next Standard®

Beyes Dental Canada Inc.
30 Casebridge Court
Toronto, Ontario, M1B 3M5
Canada
Tel: 1-416-281-1888
Email: info@beyes.ca
Website: www.beyes.ca

Your Dealer for Beyes Products



Printed in Canada.
B112EN.REV0X1716

© Beyes Canada Inc 2017. All Rights Reserved.
Specifications are subject to change without notice.
All brands are holders of their respective trademarks.