

Leading aerosol control system to protect dental professionals



DuraMax[®]

Extraoral Suction System



This aerosol control system removes contaminants produced during dental procedures including droplets, aerosols, blood, dust, and virus and uses the most advanced Quattro layer filtration system to kill virus and bacteria absorbed with over 99.9% efficiency.

During treatments, a high volume of droplets and aerosols are produced by high-speed handpieces, ultrasonic scalers, and air water syringe that can linger in the air for up to two hours. This places practitioners and office staff health at risk of infections.

DuraMax Extraoral Suction System will effectively remove these droplets and aerosols from the immediate treatment site thereby reducing the potential danger from dental aerosols and keep the air in your office safe and clean.



Remote Control

Features

EFFICIENT FILTRATION

4 layer filtration system including 1 HEPA H14 (0.1 μ m 99.99%) filtration plus germicidal UV-C

POWERFUL MOTOR

10 levels of power control with up to 270 CFM capacity

SLEEK

Compact and portable design


CONVENIENT

Unit comes with a handy remote control for smooth operation

LOW NOISE LEVEL

Built-in noise reduction technology to keep sound level below 55 decibels

Easy Filter Replacement Design

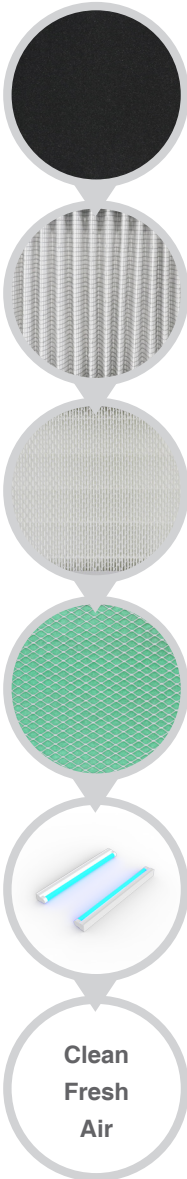

STAGE 1:
Primary Filter
Eliminate droplet, aerosol, and dust

STAGE 2:
Secondary Filter
Eliminate blood and aerosol

STAGE 3:
HEPA H14 Filter
Eliminate virus, bacteria size $\geq 0.1\mu\text{m}$ 99.99%

STAGE 4:
Activated Carbon Filter
Eliminate Odor

STAGE 5:
Germicidal Lamp UV-C
Kills viruses and bacteria



"In the context of COVID-19, airborne transmission may be possible in specific circumstances and settings in which procedures or support treatments that generate aerosols are performed"

--- World Health Organization (WHO)

Key Specifications

Model	DuraMax X1
Voltage	AC110V, 60Hz
Power	350W
Speed Options	10
Ultraviolet C (UV-C)	254nm Spectrum, 2 pieces of 8W
Maximum Motor Volume	270 CFM
Filtration	Four layers (including 1 HEPA H14 0.1µm 99.99%)
Noise Level	≤55dB
Head Diameter	2.95 inches
Arm Reach	1.4 m
Unit Weight	48 lbs
Shipping Weight	51 lbs
Shipping Dimension	L20.5 x W16 x H23.25 inches (L52 x W41 x H59 cm)
Unit Dimension	L16.5 x W12 x H20 inches (L42 x W30.5 x H46 cm)
Warranty	30 months

Order Info

Code	Model
PPS101	DuraMax X1 Extraoral Suction System, 110V
PPS102	S1 Filter
PPS103	S2 Filter
PPS104	S3 Filter
PPS105	S4 Filter
859031052	UV light tube

Under normal environmental conditions filters should be replaced every 6-12 months.

Current recommendations for filter replacement schedule for best practice may vary.



Beyes Dental Canada Inc.
 23-595 Middlefield Road
 Toronto, Ontario, M1V 3S2
 Canada
Tel: 1-416-281-1888
Email: info@beyes.ca
Website: www.beyes.ca

Your Dealer for Beyes Products

Scan the QR code for more Beyes products



Printed in Canada
 File #: ENB029
 Rev.10/09.02.21

© Beyes Canada Inc. 2020. All Rights Reserved.
 Specifications are subject to change without notice.
 All brands are holders of their respective trademarks.
 All Beyes products are manufactured in strict compliance with Health Canada regulations ensuring quality.
 Manufacturing facility is ISO 13485 certified under Canadian Medical Devices Conformity Assessment System (CMDACS).