


APEXPILOT® G3


Cordless Endo Motor

The ApexPilot G3 is a cordless and compact endodontic motor with a built-in apex locator and integrated length measuring function which makes root canal preparation safer and more efficient. The file position will be shown in an easy to read OLED screen and this endodontic motor will automatically reverse or stop as soon as the apical point is reached.

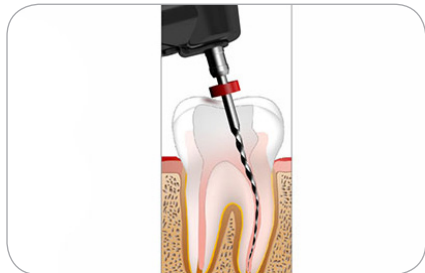

Integrated Length Determination

Automatically start when the file is inserted into the canal and stop when the file is taken out. Automatically reverse/stop when the apical point is reached.

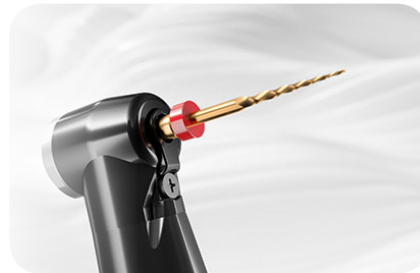

Plastic and Adjustable Contra-angle

No need to use silicone sleeve to insulate when the endo motor and apex locator are used together. Special aviation grade plastics enable contra-angle to be autoclaved to avoid cross-infection.



All-angle OLED Screen

Without viewing angle limit, clear display can be guaranteed even under strong light.



Maximum Compatibility

Wide torque and speed range make it compatible with various file system.

Key Specifications

Speed Range	100rpm~1000rpm
Torque Range	0.4-5.0N.cm
Lithium Battery	3.7V /2000mAh
Power Input	~100V-240V 50Hz/60Hz 400mA Max
Power Output	DC5V/1A
Wireless charging Frequency Range	112-205KHz
Environment Temperature	+5°C~ +40°C
Relative Humidity	30% ~ 75%
Atmospheric Pressure	70kPa ~ 106kPa
Warranty	12 Months

Order Information

Code: AL2070

Model: ApexPilot G3



Beyes Dental Canada Inc.
23-595 Middlefield Road
Toronto, Ontario, M1B 3M5
Canada
Tel: 1-855-603-1888
Email: info@beyes.ca
Web: www.beyes.ca



Lotus NL B.V.
Koningin Julianaplein 10, le Verd,
2595AA, The Hague,
Netherlands
Tel: +31645171879 (English)
+31626669008 (Dutch)
Email: info@beyes.ca

Your Dealer for Beyes Products

Scan the QR code
for more products



Printed in Canada
File #: ENB037
Rev.1/04.02.21

© Beyes Canada Inc 2021. All Rights Reserved.
Specifications are subject to change without notice.
All brands are holders of their respective trademarks.
All Beyes products are manufactured in strict compliance
with Health Canada regulations ensuring quality.
Manufacturing facility is ISO 13485 certified under Canadian
Medical Devices Conformity Assessment System (CMDCAS).